

Immediate Extraction Implantology

“ Nano-Silver Hydrogel as an Adjunct in Immediate Implant Placement: Enhancing Soft-Tissue Healing and Contour ”

Dr. Amit Gupta, MDS

Consultant Implantologist | Oral & Maxillofacial Surgeon



Profile

Dr. Amit Gupta brings extensive clinical experience in advanced implantology, including **immediate extraction implant placement, socket preservation, and esthetic soft-tissue management**. His practice emphasizes biologically guided healing and minimally invasive protocols that blend precision surgery with evidence-based regenerative materials. As an early adopter of **nano-silver hydrogel technology** in peri-implant management, Dr. Gupta has documented several successful cases demonstrating faster epithelial closure and superior soft-tissue contouring. He actively shares his insights through lectures and workshops for young clinicians, advocating a science-driven, patient-centric approach to implant dentistry.

Abstract

Immediate extraction implant placement demands precise infection control and rapid mucosal adaptation around healing abutments. Even with atraumatic extraction and ideal torque, the peri-implant mucosa often remains inflamed or irregular. In this case series, *Multident Nano-Silver Hydrogel* was used as an adjunct following immediate implant placement in two patients. The hydrogel was applied around healing abutments to contour the mucosa and promote epithelial closure. Across a **2.5-month follow-up**, both cases demonstrated excellent soft-tissue adaptation, minimal inflammation, and stable gingival collar development. These outcomes reaffirm nano-silver's dual antimicrobial and regenerative effects, translating into enhanced soft-tissue integration and improved patient comfort.





Clinical Context

Immediate implant sites carry high bacterial and inflammatory challenges. Nano-silver hydrogels provide sustained ionic release that disrupts biofilms, reduces inflammatory cytokines, and accelerates fibroblast activity. ^[1,2,3,4,5] By applying *Multident Hydrogel* directly around healing abutments, Dr. Amit Gupta achieved smoother gingival contouring and faster epithelial seal, crucial for long-term peri-implant stability.

Additionally, standard supportive therapy included painkiller tablets administered along with Multident to ensure optimal comfort and infection control.

Case 1

Immediate Extraction + Implant Placement (Maxillary Premolar & molar)


Patient Details



Age: 63y/F

Chief complaint: Cariously Fractured second molar.

Procedure

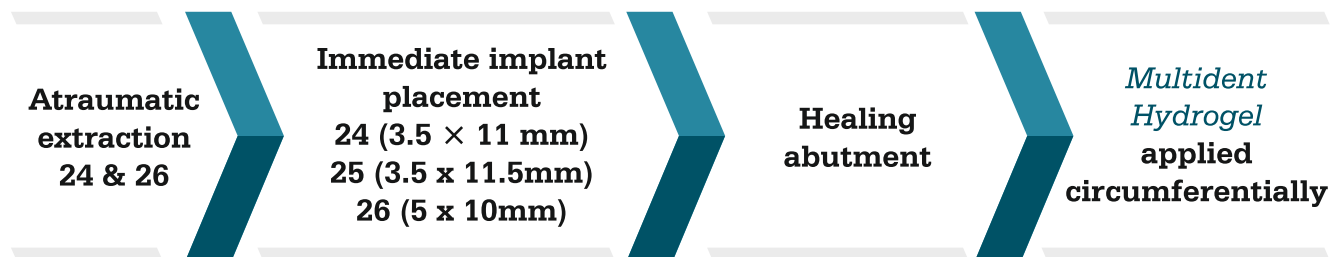

The patient was advised to use *Multident* daily for a duration of

08 weeks

Follow-up & Outcome (2.5 months)

Week 1–2

- Minimal edema
- Uneventful healing

Week 4

- Healthy attached mucosa
- Stable epithelial maturation
- Light-pink tone

Week 8–10 (2.5 months)

- Fully mature gingival collar with ideal emergence profile
- No mucosal tenderness or inflammation
- Provisionalization feasible with highly stable peri-implant tissues



Clinical Insight

Nano-silver reduced microbial irritation at the mucosal margin and supported soft-tissue sculpting.

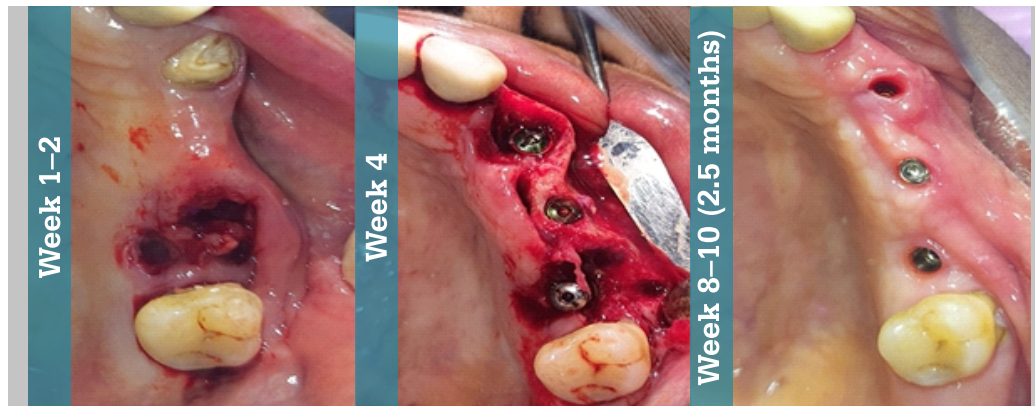


Figure 1: Immediate extraction and implant placement in the maxillary premolar region.

Case 2

Immediate Extraction + Implant Placement (Mandibular Molar)


Patient Details



Age: 54y/F

Chief complaint: Fractured first premolar and grade 1 mobility with first molar.

Procedure

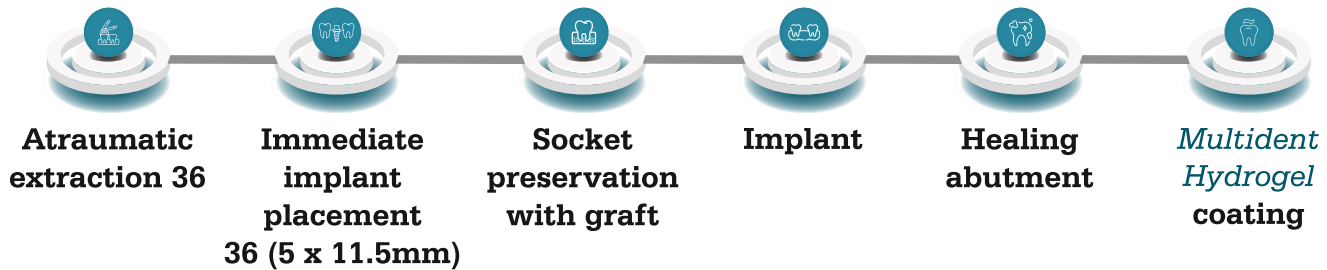

The patient was advised to use Multident daily for a duration of

08
weeks

Follow-up & Outcome (2.5 months)

Week 1-2

- Minimal discomfort
- Early fibrin integrity intact

Week 4

- Near-complete epithelial coverage
- No peri-implant erythema

Week 8-10 (2.5 months)

- Fully keratinized mucosal seal
- No bleeding on probing
- **Radiograph:** clear crestal bone line and no soft-tissue gap





Clinical Insight

In grafted sites, nano-silver promoted faster mucosal sealing and easier prosthetic impression timing.



Figure 2: Immediate extraction with socket grafting and implant placement in the mandibular molar region.

Discussion & Expert Commentary

Silver nanoparticles have documented antibacterial activity against *S. mutans*, *P. gingivalis*, and *F. nucleatum*,^[1,2,3] they also modulate inflammatory mediators (IL-1 β , TNF- α).^[4,5] Their collagen-stimulatory effect aids epithelial re-attachment and contour formation around abutments. In both cases, Dr. Gupta noted faster healing, reduced bleeding, and more predictable esthetics-outcomes consistent with published nano-silver evidence.

Mechanism Summary

Mechanism	Molecular Effect	Clinical Outcomes
Antimicrobial	Destroys peri-implant biofilm ^[1,2,3]	Reduces infection risk
Anti-inflammatory	Suppresses IL-1 β , TNF- α ^[4,5]	Minimizes edema & soreness
Regenerative	Stimulates fibroblast & collagen synthesis ^[5]	Faster epithelial seal

Key Takeaways

1. *Multident Nano-Silver Hydrogel* enhances soft-tissue healing after immediate implant placement.
2. Offers antimicrobial protection without tissue irritation.
3. Improves gingival collar formation and peri-implant esthetics within a clinically relevant prosthetic phase (8–10 weeks).



Multident™



References: 1) Raghav P, Khara A K, Bisht S. Comparative evaluation of antimicrobial properties of silver nanoparticles and chlorhexidine mouthwashes. *J Clin Orthod Dent.* 2023. 2) Nguyen T et al. Comparing the effectiveness of 0.12% chlorhexidine and silver nanoparticle mouthwashes. *Vietnam J Dent Sci.* 2024. 3) Rusu A et al. Nanoparticles in periodontitis therapy: A review of current trends. *Int J Mol Sci.* 2023. 4) Hernández-Sierra JF et al. Bactericidal activity of silver nanoparticles on oral biofilms related to periodontal disease. *Bioengineering.* 2022; 14(6):311. 5) Chen S et al. Comparison of antimicrobial and wound-healing effects of silver nanoparticle and chlorhexidine mouthwashes. *Clin Oral Invest.* 2022.

Disclaimer: The case reports presented herein are intended for educational purposes only and reflect the clinical opinions and experiences of the contributing experts. They do not constitute medical advice, endorsement, or a standard of care. Healthcare professionals must use their own clinical judgment when applying this information to patient care. Individual patient circumstances may vary, and consultation with appropriate specialists is always recommended.

